

## PROTOCOL

### Murine IFN- $\gamma$ / IL-5 Double-Color FluoroSpot Assay

#### CONTENTS

##### **Murine IFN- $\gamma$ /IL-5 Capture Kit:**

- Murine IFN- $\gamma$  Capture Ab
- Murine IL-5 Capture Ab
- CTL-Test™ Medium
- Diluent A
- Diluent B
- Diluent C
- Plates: 96-well,  
low auto-fluorescent,  
high-protein-binding,  
PVDF filter plates
- Adhesive plate  
sealing sheet
- Protocol

##### **Murine IFN- $\gamma$ Detection Kit:**

- Anti-murine IFN- $\gamma$   
(FITC) Detection Ab
- Anti-FITC Alexa Fluor®  
488

##### **Murine IL-5 Detection Kit:**

- Anti-murine IL-5  
(Biotin) Detection Ab
- Strep CTL-Red™



Visit our YouTube channel  
for several helpful videos on  
working with ELISPOT and  
FluoroSpot assays:  
[www.youtube.com/user/  
ImmunoSpot](http://www.youtube.com/user/ImmunoSpot).

Protocols and technical  
resources available at  
[www.immunospot.com](http://www.immunospot.com).

#### PROCEDURE

##### DAY 0 — STERILE CONDITIONS

- Prepare *Murine IFN- $\gamma$ /IL-5 Capture Solution* and prepare 70% ethanol (see Solutions).
- Remove plate underdrain, pipette 150 $\mu$ l of 70% ethanol into each well and incubate for less than one minute. Add 150 $\mu$ l of PBS, decant, and wash with 150 $\mu$ l of PBS two more times. (If using strip plates, there is no underdrain to remove before prewetting.) **Note:** Activation of the membrane with ethanol is instantaneous and can be seen visually as a graying of the membrane. Ethanol should be washed off as quickly as possible following activation.
- Replace underdrain and immediately (before plate dries) pipette 80 $\mu$ l/well *Murine IFN- $\gamma$ /IL-5 Capture Solution*. Seal plate with parafilm and incubate at 4°C overnight.

##### DAY 1 — STERILE CONDITIONS


- Prepare *CTL-Test™ Medium* (see Solutions).
- Prepare antigen/mitogen solutions at two times final concentration in *CTL-Test™ Medium*.
- Decant plate containing *Capture Solution* from Day 0 and wash one time with 150 $\mu$ l PBS.
- Plate antigen/mitogen solutions, 100 $\mu$ l/well. Ensure the pH and temperature are ideal for cells by placing the plate containing antigens into a 37°C incubator for 10-20 minutes before plating cells.
- Adjust cells to desired concentration in *CTL-Test™ Medium*, e.g.: 3 million/ml corresponding to 300,000 cells/well (cell numbers can be adjusted according to expected spot counts since 100,000-800,000 cells/well will provide linear results). Keep cells at 37°C in humidified incubator, 9% CO<sub>2</sub> while processing cells and until plating.
- Plate cells 100 $\mu$ l/well using large orifice tips. Once completed, gently tap the sides of the plate and immediately place into a 37°C humidified incubator, 9% CO<sub>2</sub>.
- Incubate for 24 hours. Do not stack plates. Avoid shaking plates by carefully opening and closing incubator door. Do not touch plates during incubation.

##### DAY 2

- Prepare Wash Buffers for the day: PBS, distilled water and Tween-PBS.
- Prepare *Anti-murine IFN- $\gamma$ /IL-5 Detection Solution* (see Solutions).
- Wash plate two times with PBS and then two times with 0.05% Tween-PBS, 200 $\mu$ l/well each time.
- Add 80 $\mu$ l/well *Anti-murine IFN- $\gamma$ /IL-5 Detection Solution*. Incubate at room temperature, two hours.
- Prepare *Tertiary Solution* (see Solutions).
- Wash plate three times with 0.05% Tween-PBS, 200 $\mu$ l/well.
- Add 80 $\mu$ l/well of *Tertiary Solution*. Incubate at room temperature, one hour.
- Decant and wash plate three times with distilled water, 200 $\mu$ l/well. (Optimal results can be seen when the last water wash is filtered through with a vacuum manifold to get rid of any unbound tertiary.)
- Remove protective underdrain from the plate and rinse back of plate with tap water.
- Air-dry plate for two hours in running laminar flow hood or for 24 hours face down on paper towels on bench top.
- Scan and count plate. (CTL has scanning and analysis services available and offers a trial version of ImmunoSpot® Software with the purchase of any kit. Email [kitscanningservices@immunospot.com](mailto:kitscanningservices@immunospot.com).) **Note:** Fluorescent signals must be read with compatible light source(s) and filter sets. The optimized settings differ depending on the model of instrument used. Please consult with Technical Support for assistance at +1 216-791-5084.

See other side for Solutions and Technical Tips.

For laboratory research use only. Not for use in diagnostic or therapeutic procedures.

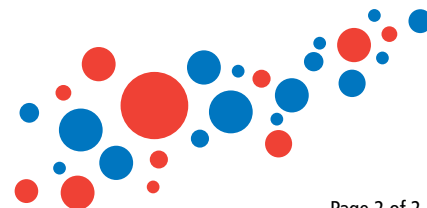

## SOLUTIONS

**All solutions should be freshly-made prior to use. It is important to quick-spin the vials before use to ensure content volumes.**

- **70% Ethanol (not included):** Dilute 190-200 proof ethanol. For 10ml, add 7ml of ethanol to 3ml of distilled water.
- **CTL-Test™ Medium:** Prepare medium by adding 1% fresh L-glutamine. Amount of medium needed will depend on variables such as cell yield and number of samples tested but will be no less than 20ml for one full plate.
- **Capture Solution:** Dilute *Murine IFN-γ* and *Murine IL-5 Capture Antibodies* in *Diluent A*. For one plate, add 60μl of *Murine IFN-γ* and 80μl of *Murine IL-5 Capture Antibodies* to 10ml of *Diluent A*.
- **Detection Solution:** Dilute *Anti-murine IFN-γ (FITC)* and *Anti-murine IL-5 (Biotin) Detection Antibodies* in *Diluent B* and filter. For one plate, add 20μl of *Anti-murine IFN-γ (FITC)* and 40μl of *Anti-murine IL-5 (Biotin) Detection Antibodies* to 10ml of *Diluent B* and filter through a 0.22μm filter.
- **Tertiary Solution:** Dilute *Anti-FITC Alexa Fluor® 488* (visualizes IFN-γ) and *Strep CTL-Red™* (visualizes IL-5) and in *Diluent C* filter. For one plate, add 25μl of *Anti-FITC Alexa Fluor® 488* and 50μl of *Strep CTL-Red™* in 10ml of *Diluent C* and filter through a 0.22μm filter.

### Wash Buffers (not included)

For each plate prepare:

- 0.05% Tween-PBS: 100μl Tween-20 in 200ml PBS
- PBS, sterile, 100ml
- Distilled water, 100ml

## TECHNICAL TIPS

- Upon successful completion of the assay and proper excitation, IFN-γ spots will fluoresce green, IL-5 spots will fluoresce red.
- To maximize the use of each plate, an adhesive plate-sealing sheet has been included that can be adhered to the top of the plate to cover unused wells for use in subsequent assays. Use your thumbs to firmly adhere the sheet to the plate and a razor blade to cut the sheet to expose only the necessary wells.
- We highly recommend the use of CTL Serum-free Media for freezing, washing, and testing cells. Even brief exposure to a mitogenic serum can cause high background while other sera can have suppressive effects. CTL also recommends using the CTL-LDC™ Kit for accurate live/dead cell counts.
- Deviations from specified temperatures, timing requirements, number of washing steps, and specified reagent preparation volumes may alter the performance of the assay.
- Plates may be washed manually or with a suitable automated plate washer with adjusted pin length and flow rate so membranes and spots are not damaged (CTL recommends the CTL 405LSR).
- To avoid damage to the PVDF membrane in the wells, do not touch the membrane with pipette tips or with the plate washer. The PVDF membrane is permeable and protected by an underdrain. Avoid direct contact between the well bottom and wet surfaces, including paper towels or any other materials that will absorb liquid.
- While processing plates, the PVDF membrane at the bottom of the wells must remain wet.
- When underdrain and gloves are wet, the underdrain may be slippery and difficult to remove. Wipe gloves and underdrain with paper towel before removing.
- After completion of the experiment, do not dry the ELISPOT assay plates at temperatures exceeding 37°C as this may cause the membrane to crack.
- Spots may not be readily visible while the membrane is still wet. Scan and count plates only after membranes have completely dried.
- Higher background appearing in the control wells can be potentially overcome using the following steps:
  - When working with precultured cells, wash the cells thoroughly in CTL-Wash™ prior to the experiment in order to avoid carryover of cytokines and other substances; use CTL-Test™ for testing cells.
  - The SmartCount™ module of the ImmunoSpot® counting software automatically recognizes spots over high background or uneven background, correcting background deviations. The Autogating™ module will help discern between T cell-derived and background spots. The CTL technical support team will gladly assist you with using the ImmunoSpot® Software for the analysis of complicated test results.
- Data analysis: The CTL ImmunoSpot® Analyzers along with the ImmunoSpot® Software have advanced features that permit automated, objective recognition of spots, gating and counting. An ELISPOT data management tool, SpotMap®, is also available to facilitate high-throughput ELISPOT work.

The CTL team will gladly assist you with data analysis and troubleshooting, as well as in customizing ELISPOT assays to suit your needs. Please contact us at [kits@immunospot.com](mailto:kits@immunospot.com).



Visit our YouTube channel for several helpful videos on working with ELISPOT and FluoroSpot: [www.youtube.com/user/ImmunoSpot](http://www.youtube.com/user/ImmunoSpot).

Protocols and technical resources are available at [www.immunospot.com](http://www.immunospot.com).

See other side for Contents and Procedure.

For laboratory research use only. Not for use in diagnostic or therapeutic procedures.



**Cellular Technology Limited** • 20521 Chagrin Blvd. • Shaker Hts., OH 44122-5350 USA  
+1 216-791-5084 Phone • +1 888-791-4005 Toll Free US • +1 216-751-1373 Fax  
[kits@immunospot.com](mailto:kits@immunospot.com) • [www.immunospot.com](http://www.immunospot.com)

**CTL-Europe GmbH** • Hans-Boeckler-Str. 19-29, 53225 • Bonn, Germany  
+49 (0) 228 8544 78-0 Phone • +49 (0) 228 8544 78-22 Fax  
[ctl-europe@immunospot.eu](mailto:ctl-europe@immunospot.eu) • [www.immunospot.eu](http://www.immunospot.eu)